



**OFFICE OF STATEWIDE HEALTH PLANNING AND DEVELOPMENT  
FACILITIES DEVELOPMENT DIVISION**

**APPLICATION FOR OSHPD PREAPPROVAL  
OF MANUFACTURER'S CERTIFICATION (OPM)**

OFFICE USE ONLY

APPLICATION #: OPM-0207-13

**OSHPD Preapproval of Manufacturer's Certification (OPM)**

Type:  New  Renewal  Update to Pre-CBC 2013 OPA Number: \_\_\_\_\_

**Manufacturer Information**

Manufacturer: CareFusion

Manufacturer's Technical Representative: Jared Zamaloff

Mailing Address: 10020 Pacific Mesa Blvd., San Diego, CA. 92121

Telephone: (818) 876-4287

Email: Jared.zamaloff@carefusion.com

**Product Information**

Product Name: Single Storage Cabinet

Product Type: Other Electrical and Mechanical Components

Product Model Number: Pyxis SupplyStation, Pyxis MedStation and Pyxis Medstation Console

General Description: The Pyxis Supply Station – system provides easy access to needed supplies on nursing floors and throughout your healthcare facility. This secure storage device provides your staff with the ability to document supply usage, in real-time. Pyxis MedStation system is a leading automated dispensing system supporting decentralized medication management. Pyxis MedStation Console-server collects data from Pyxis supply management systems located throughout your facility.

**Applicant Information**

Applicant Company Name: EASE Co.

Contact Person: Jonathan Roberson, S.E.

Mailing Address: 5877 Pine Ave. Suite 210, Chino Hills, CA. 91709

Telephone: (909) 606-7622

Email: J.Roberson@EASECo.com

I hereby agree to reimburse the Office of Statewide Health Planning and Development review fees in accordance with the California Administrative Code, 2013.

Signature of Applicant: \_\_\_\_\_

Date: 3/30/15

Title: Principal Engineer

Company Name: EASE Co.

"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs"

STATE OF CALIFORNIA – HEALTH AND HUMAN SERVICES AGENCY  
OSH-FD-700 (REV 1/24/13)



osHPD

Page 1 of 2



**OFFICE OF STATEWIDE HEALTH PLANNING AND DEVELOPMENT  
FACILITIES DEVELOPMENT DIVISION**

**Registered Design Professional Preparing Engineering Recommendations**

Company Name: EASE Co.

Name: Jonathan Roberson, S.E. California License Number: S4197

Mailing Address: 5877 Pine Ave. Suite 210, Chino Hills, CA. 91709

Telephone: 909-606-7667 Email: J.Roberson@EASECo.com

**OSHDP Special Seismic Certification Preapproval (OSP)**

Special Seismic Certification is preapproved under OSP-  
(Separate application for OSP is required)

Special Seismic Certification is not preapproved

**Certification Method(s)**

Testing in accordance with:  ICC-ES AC156  FM 1950-10

Other\* (Please Specify): \_\_\_\_\_

\*Use of test criteria other than those adopted by the California Building Standards Code, 2013 (CBSC 2013) for component supports and attachments are not permitted. For distribution system, interior partition wall, and suspended ceiling seismic bracings, test criteria other than those adopted in the CBSC 2013 may be used when approved by OSHDP prior to testing.

Analysis

Experience Data

Combination of Testing, Analysis, and/or Experience Data (Please Specify): \_\_\_\_\_

**List of Attachments Supporting the Manufacturer's Certification**

Test Report  Drawings  Calculations  Manufacturer's Catalog

Other(s) (Please Specify): \_\_\_\_\_

**OFFICE USE ONLY – OSHDP APPROVAL VALID FOR CBC 2013 ONLY**

Signature:  Date: 06/19/2015

Print Name: Jeffrey Kikumoto

Title: SSE

Condition of Approval (if applicable): \_\_\_\_\_

"Access to Safe, Quality Healthcare Environments that Meet California's Diverse and Dynamic Needs"





**EQUIPMENT ANCHORAGE  
& SEISMIC ENGINEERING**

5877 Pine Ave, Ste. 210  
Chino Hills, CA. 91709  
Phn: (909) 606-7622

Office of Statewide Health Planning and Development  
**PREAPPROVAL OF MANUFACTURER'S CERTIFICATION**  
**OPM-0207-13**

**THIS PREAPPROVAL CONFORMS TO THE 2013 CALIFORNIA BUILDING CODE**

MANUFACTURER: **CAREFUSION**  
EQUIPMENT NAME: **SINGLE STORAGE CABINETS**

Sheet: 1 of 12  
Date: 6/18/15

**GENERAL NOTES**

1. THIS OSHPD PREAPPROVAL OF MANUFACTURER'S CERTIFICATION (OPM) IS BASED ON THE 2013 CBC. THE DEMANDS (DESIGN FORCES) FOR USE WITH THIS OPM SHALL BE BASED ON THE 2013 CBC
2. THIS DOCUMENT MAY ONLY BE USED WITH THE EXPRESS WRITTEN CONSENT OF THE MANUFACTURER LISTED ABOVE FOR THE SPECIFIC PROJECT SITE AND INSTALLATION LOCATION. THIS DOCUMENT IS INVALID WITHOUT SUCH CONSENT.
3. THIS PREAPPROVAL CONFORMS TO THE 2013 CALIFORNIA BUILDING CODE WHERE  $S_{ds}$  IS NOT GREATER THAN 1.15 & 2.20. SEE DETAILS FOR APPLICABILITY.
4. FORCES PER ASCE 7-10 SECTION 13.3.1, EQUATIONS 13.3-1, 13.3-2 & 13.3-3,  
WHERE  $S_{ds} = 1.15$ ,  $a_p = 1.0$ ,  $I_p = 1.5$ ,  $R_p = 1.5$ ,  $z/h = 0$  AT CONCRETE SLAB. SEE FOLLOWING SHEETS FOR  $\Omega$ .  
WHERE  $S_{ds} = 2.20$ ,  $a_p = 1.0$ ,  $I_p = 1.5$ ,  $R_p = 1.5$ ,  $z/h \leq 1$  AT CONCRETE SLAB ON METAL DECK. SEE FOLLOWING SHEETS FOR  $\Omega$ .
5. THIS PREAPPROVAL COVERS ONLY THE SUPPORTS AND ATTACHMENTS OF THE EQUIPMENT TO THE STRUCTURE.
6. ALL DESIGN FORCES SHOWN ON THE DRAWINGS ARE FACTORED LOADS THAT SHALL BE USED FOR STRENGTH DESIGN.
7. CONCRETE SLAB ON METAL DECK DETAIL VALID FOR DEMANDS SHOWN AT ANY ELEVATION IN THE BUILDING. (i.e.  $z/h \leq 1$ )
8. CONCRETE SLAB ON GRADE DETAIL VALID FOR DEMANDS SHOWN AT ANY ELEVATION BELOW GRADE. (i.e.  $z/h = 0$ )
9. **RESPONSIBILITIES OF THE STRUCTURAL ENGINEER OF RECORD OF THE BUILDING**

- A. PROVIDE SUPPORTING STRUCTURE TO SUPPORT WEIGHTS AND FORCES SHOWN IN ADDITION TO ALL OTHER LOADS.
- B. VERIFY THAT THE INSTALLATION IS IN CONFORMANCE WITH THE 2013 CBC AND WITH THE DETAILS, MATERIAL AND GAGE OF THE UNIT WHERE ATTACHMENTS ARE MADE AGREE WITH THE INFORMATION SHOWN ON THE PREAPPROVAL DOCUMENTS.
- C. VERIFY THAT PROJECT SPECIFIC VALUES OF  $S_{ds}$  &  $z/h$  RESULT IN SEISMIC FORCES ( $E_h$ ,  $E_v$ ) THAT DO NOT EXCEED THE VALUES ON THE DETAILS.
- D. VERIFY THAT THE CONCRETE SLAB TO WHICH THE EQUIPMENT IS ANCHORED MEETS THE REQUIREMENTS OF THE APPLICABLE ICC ESR.
- E. VERIFY THAT THE ANCHORS ARE AN ADEQUATE DISTANCE FROM ANY SLAB EDGES OR OPENINGS (SEE TYPICAL DETAIL ON SHEET 2).
- F. VERIFY THAT ALL NEW OR EXISTING ANCHORS ARE AN ADEQUATE DISTANCE FROM THE UNIT ATTACHMENTS AND CHECK FOR INTERACTION WHERE OTHER ANCHORS ARE WITHIN 18" OR  $6h_{ef}$  FROM THIS UNIT'S ANCHORS.



## CAREFUSION

### SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**

SHEET

**2**

OF **12** SHEETS

#### 10. EXPANSION ANCHORS:

A. ATTACHMENT IS TO BE MADE WITH THE ANCHORS LISTED BELOW AND INSTALLED AS DESCRIBED IN THE CORRESPONDING ICC REPORT.

Anchor Diameter	Concrete Type	Min. f'c (psi)	Anchor Type	ICC Report No.	Min. Embed.	Min. Spacing	Min. Edge Dist.	Min. Conc. Thickness	Torque Test	Direct Tension
3/8"	Sand Light Weight	3000	Hilti Kwik Bolt TZ	ESR-1917	2"	N/A	N/A	See Sheet 10 of 12	25 FT-LB	1186 lb
3/8"	Normal Weight	3000	Hilti Kwik Bolt TZ	ESR-1917	2"	8"	12"	4"	25 FT-LB	1515 lb

B. THIS PREAPPROVAL ALLOWS FOR UP TO A MAXIMUM OF 2 ADJACENT CONCRETE SLAB EDGES, 12" AWAY MINIMUM (i.e. - CORNER). SEE ADJACENT DETAIL FOR ADDITIONAL MINIMUM ALLOWABLE CONCRETE EDGE DISTANCES.


C. TESTING OF EXPANSION ANCHORS PER 2013 CBC, 1913A.7: TESTING SHALL BE DONE IN THE PRESENCE OF THE SPECIAL INSPECTOR AND A REPORT OF THE TEST RESULTS SHALL BE SUBMITTED TO OSHPD

(i) AFTER AT LEAST 24 HOURS HAVE ELAPSED SINCE INSTALLATION, DIRECT PULL TENSION TEST OR TORQUE TEST AT LEAST 50% OF THE ANCHORS.

(ii) ACCEPTANCE CRITERIA:

- DIRECT TENSION TEST: THE ANCHOR SHOULD HAVE NO OBSERVABLE MOVEMENT AT THE TEST LOAD. A PRACTICAL WAY TO DETERMINE OBSERVABLE MOVEMENT IS THAT THE WASHER BECOMES LOOSE.
- TORQUE TEST: THE APPLICABLE TORQUE MUST BE ACHIEVED WITHIN THE FOLLOWING LIMITS: WEDGE TYPE : 1/2 TURN OF THE NUT

(iii) IF ANY ANCHOR FAILS, TEST ALL ANCHORS.

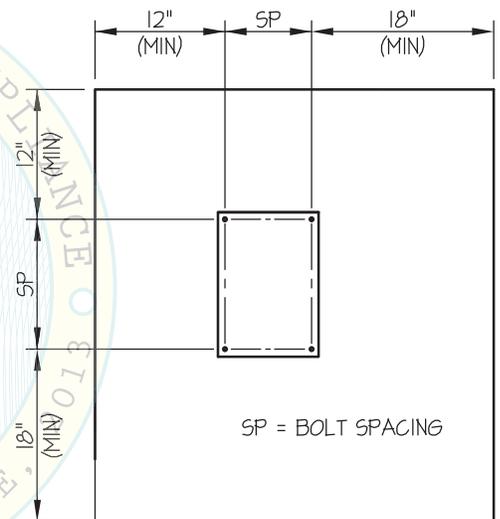

TYPICAL CONCRETE EDGE DETAIL

#### 11. BOLTS THROUGH CONCRETE ON METAL DECK

- A. BOLTS SHALL BE TORQUED BY 3/4 TURN OF THE NUTS AFTER THE SNUG TIGHT (THE SNUG-TIGHT CONDITION IS DEFINED AS THE TIGHTNESS REQUIRED TO BRING THE CONNECTED PLIES INTO FIRM CONTACT) CONDITION IS ACHIEVED, UNLESS OTHERWISE NOTED.
- B. THROUGH BOLT HOLES SHALL BE 1/16" LARGER THAN BOLT SIZE (HOLE SIZE = BOLT SIZE + 1/16) FOR CONCRETE.
- C. THROUGH-BOLTS IN CONCRETE SHALL RECEIVE SPECIAL INSPECTION AND TESTING (THROUGH BOLTS WITH STEEL TO STEEL CONNECTION IN TENSION DO NOT REQUIRE TENSION TESTING) IN ACCORDANCE WITH REQUIREMENTS FOR POST-INSTALLED ANCHORS.



## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**


SHEET

**3**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE SLAB

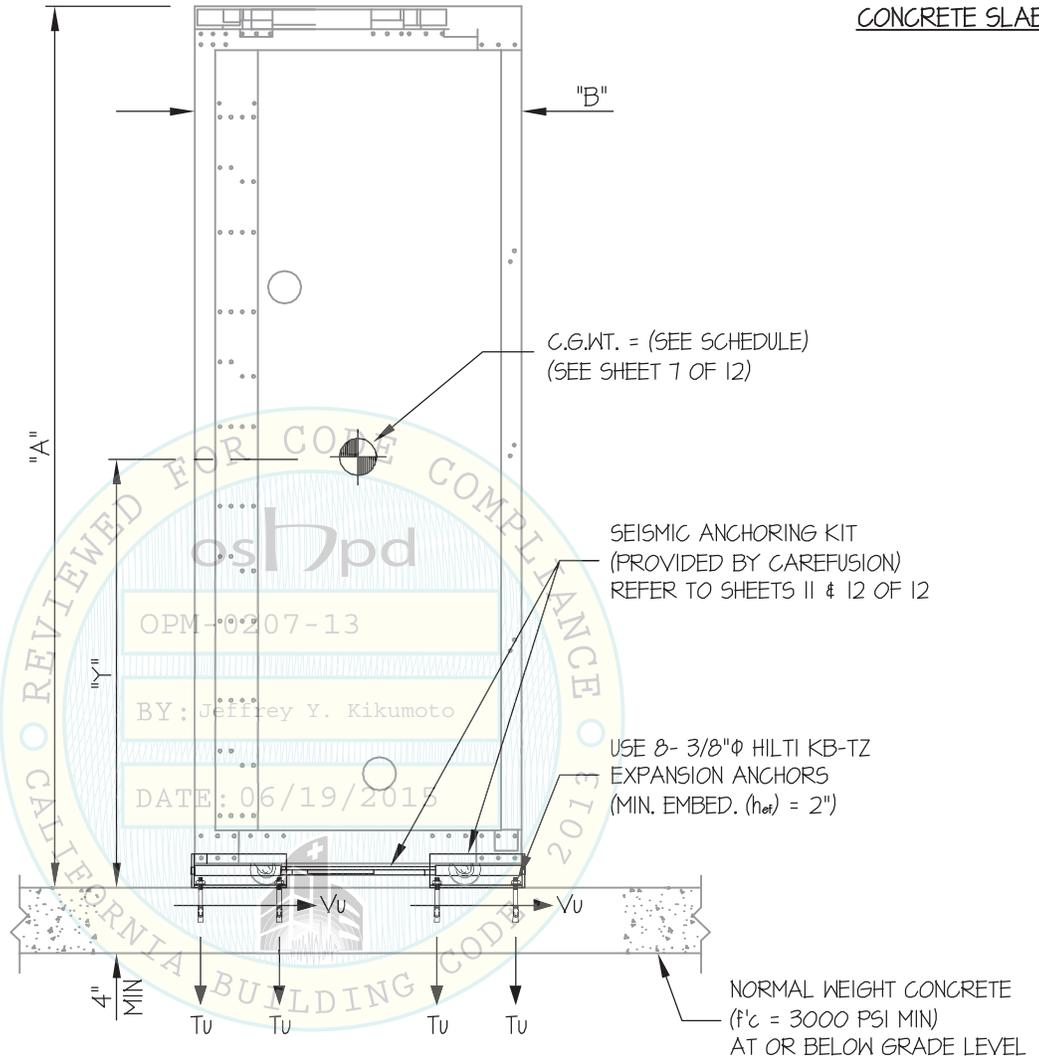

**FRONT ELEVATION**

SEE SHEET 6 OF 12 FOR  
WEIGHTS, DIMENSIONS, & BOLT FORCES

NOTES:

- FORCES ARE DETERMINED PER 2013 CALIFORNIA BUILDING CODE AND ASCE 7-10. STRENGTH DESIGN IS USED. ( $\alpha_p = 1.0$ ,  $l_p = 1.5$ ,  $R_p = 1.5$ ,  $\Omega_0 = 1.5$ ,  $z/h = 0$ )
- CENTER OF GRAVITY (C.G.) AND WEIGHT ARE THE GOVERNING PARAMETERS FOR DESIGN. THIS PREAPPROVAL ENCOMPASSES ALL WEIGHTS UP TO THE MAXIMUM WEIGHT SHOWN.
- STRUCTURAL ENGINEER OF RECORD FOR THE BUILDING SHALL PROVIDE SUPPORT STRUCTURE DESIGNED TO SUPPORT WEIGHTS AND FORCES SHOWN IN COMBINATION WITH ALL OTHER LOADS THAT MAY BE PRESENT.
- SEE GENERAL NOTES: SHEETS 1 AND 2.



## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**

SHEET

**4**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

MAX  $S_{Ds} \leq 1.15$

CONCRETE SLAB

USE 8- 3/8"  $\phi$  HILTI KB-TZ  
EXPANSION ANCHORS  
(MIN. EMBED. ( $h_{dev}$ ) = 2")

NOTE: SEE SHEET 12 OF 12  
FOR G1 & G2 DIMENSIONS



C.G. WT. = (SEE SCHEDULE)  
( $\bar{Y}$  = "Y")  
(SEE SHEET 6 OF 12)

STEEL DOWEL  
(1/2" DIA, ASTM A1018 HSLAS,  
CLASS I, GRADE 50,  $F_y=50$  KSI (MIN))  
(EA SIDE)

DOWEL EXTEND A MIN 1/8"  
PAST FACE OF BRACKET  
(TYP)

RIGHT & LEFT  
FLOOR BRACKETS  
(1/4" THK, ASTM A1018 HSLAS,  
CLASS I, GRADE 50,  $F_y=50$  KSI (MIN))

**PLAN AT BASE**  
(SYSTEM EC SINGLE UNIT SHOWN)

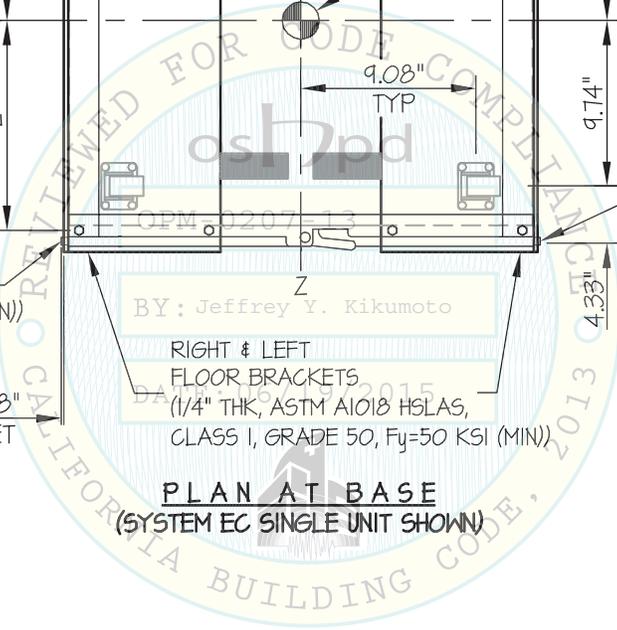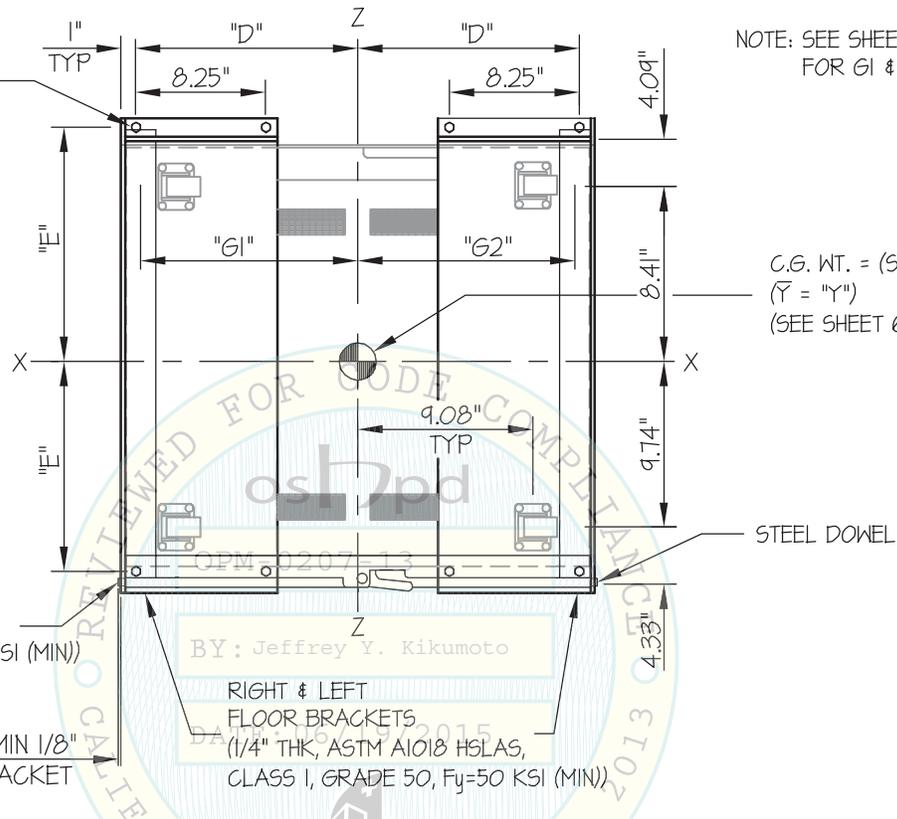

## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**

SHEET

**5**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

1.15 < MAX S<sub>D5</sub> ≤ 2.20

CONCRETE SLAB

USE 8- 3/8"Ø HILTI KB-TZ  
EXPANSION ANCHORS  
(MIN. EMBED. (h<sub>ef</sub>) = 2")

NOTE: SEE SHEET 12 OF 12  
FOR G1 & G2 DIMENSIONS

C.G. WT. = (SEE SCHEDULE)  
(Y = "Y")  
(SEE SHEET 6 OF 12)



STEEL DOWEL  
(1/2" DIA, ASTM A1018 HSLAS,  
CLASS I, GRADE 50, F<sub>y</sub>=50 KSI (MIN))  
(EA SIDE)

DOWEL EXTEND A MIN 1/8"  
PAST FACE OF BRACKET  
(TYP)

BY: Jeffrey Y. Kikumoto

RIGHT & LEFT  
FLOOR BRACKETS  
(1/4" THK, ASTM A1018 HSLAS,  
CLASS I, GRADE 50, F<sub>y</sub>=50 KSI (MIN))

**PLAN AT BASE**  
(CONSOLE UNIT SHOWN)

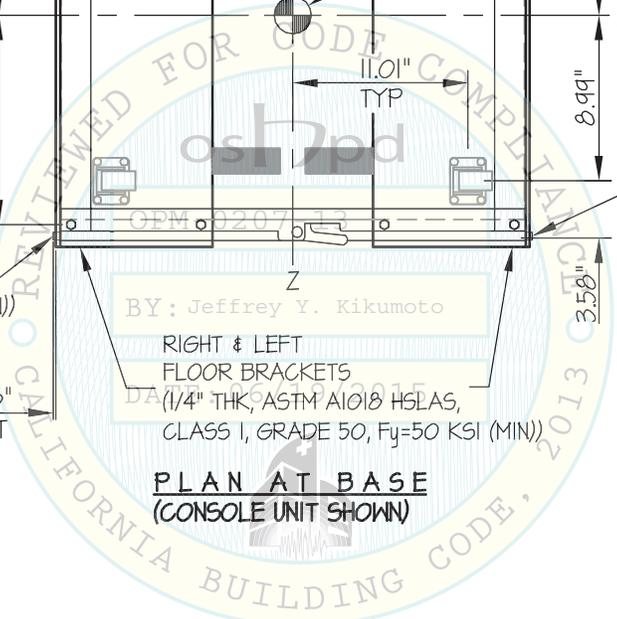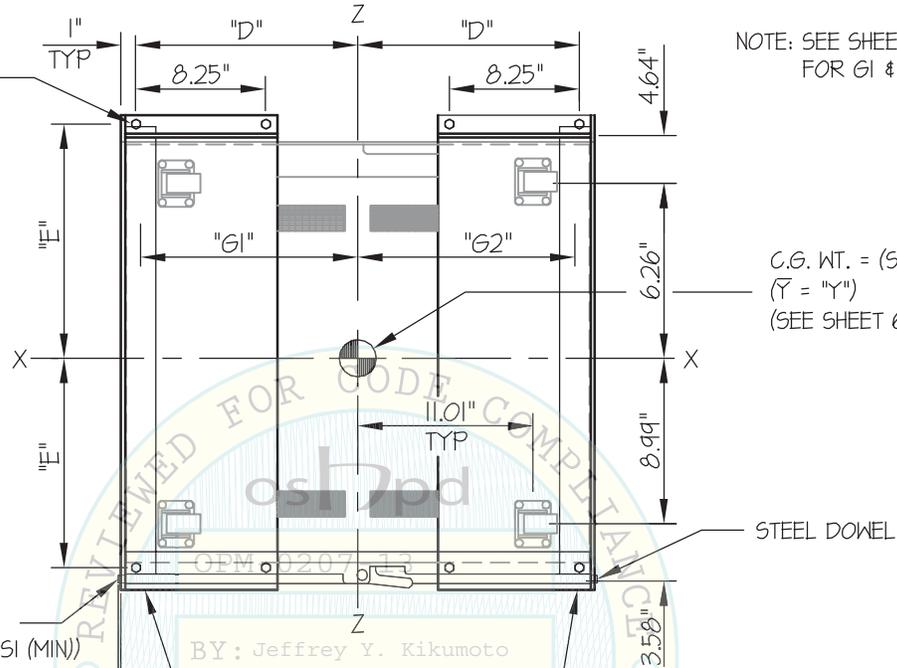

*Jonathan Roberson*  
REGISTERED PROFESSIONAL ENGINEER  
No. 4197  
EXP. 6-30-2016  
6/18/15  
STATE OF CALIFORNIA

## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**

SHEET

**6**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE SLAB


MAX Sds	MODEL	WEIGHT (LB)	"A" (in)	"Y" (in)	"B" (in)	"D" (in)	"E" (in)	* TU (lb)	* VU (lb)
2.20	CONSOLE UNIT	360	48.50	24.25	30.00	14.375	11.75	420	124
1.15	MED STATION	1050	46.75	23.38	22.75	10.75	13.25	565	214
1.15	PYXIS STATION SINGLE UNIT	1460	78.06	39.03	30.30	14.53	13.25	1138	268
1.15	SYSTEM 30 SINGLE UNIT	1460	78.06	39.03	30.30	14.53	13.25	1138	268
1.15	SYSTEM EC SINGLE UNIT	1460	78.06	39.03	30.30	14.53	13.25	1138	268

\* VALUES INCLUDE  $\Omega$

NOTE: VALUES ARE FOR  $z/h = 0$ .

BY: Jeffrey Y. Kikumoto

DATE: 06/19/2015

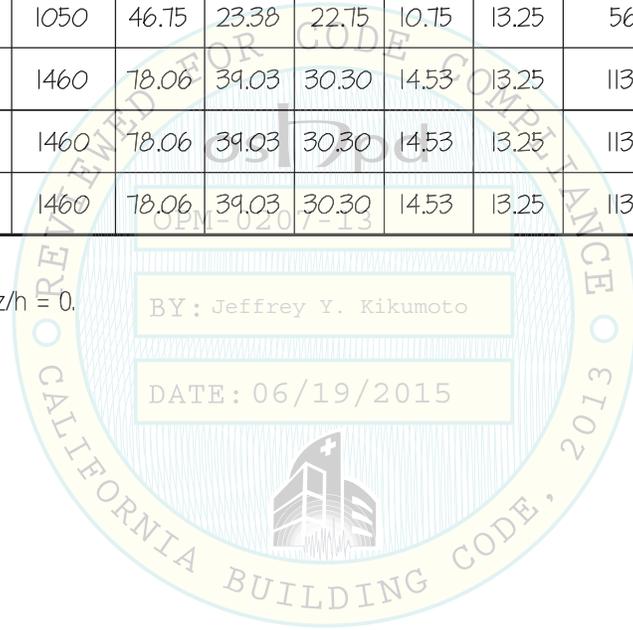

## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**


SHEET

**7**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE SLAB ON METAL DECK

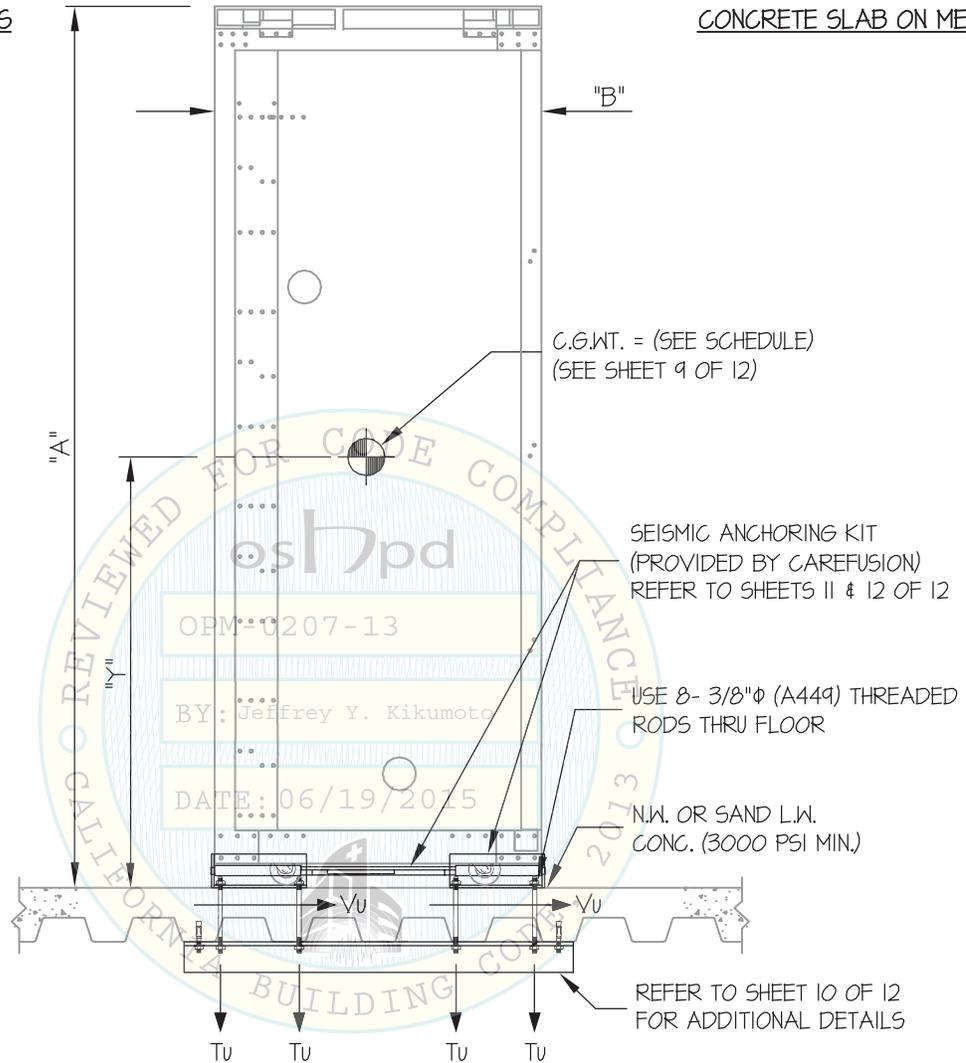

**FRONT ELEVATION**

SEE SHEET 9 OF 12 FOR WEIGHTS, DIMENSIONS, & BOLT FORCES

**NOTES:**

- FORCES ARE DETERMINED PER 2013 CALIFORNIA BUILDING CODE AND ASCE 7-10. STRENGTH DESIGN IS USED. ( $S_{bs} = 2.20$ ,  $a_p = 1.0$ ,  $I_p = 1.5$ ,  $R_p = 1.5$ ,  $\Omega_o = 1.5$ ,  $z/h \leq 1$ )

HORIZONTAL FORCE ( $E_h$ ) =  $2.64 W_p$

HORIZONTAL FORCE ( $E_{mh}$ ) =  $3.96 W_p$  (FOR CONCRETE ANCHORAGE)

VERTICAL FORCE ( $E_v$ ) =  $0.44 W_p$

- CENTER OF GRAVITY (C.G.) AND WEIGHT ARE THE GOVERNING PARAMETERS FOR DESIGN. THESE CALCULATIONS ENCOMPASS ALL WEIGHTS UP TO THE MAXIMUM WEIGHT SHOWN.
- STRUCTURAL ENGINEER OF RECORD FOR THE BUILDING SHALL PROVIDE SUPPORT STRUCTURE DESIGNED TO SUPPORT WEIGHTS AND FORCES SHOWN IN COMBINATION WITH ALL OTHER LOADS THAT MAY BE PRESENT.



## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**


SHEET

**8**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE SLAB ON METAL DECK

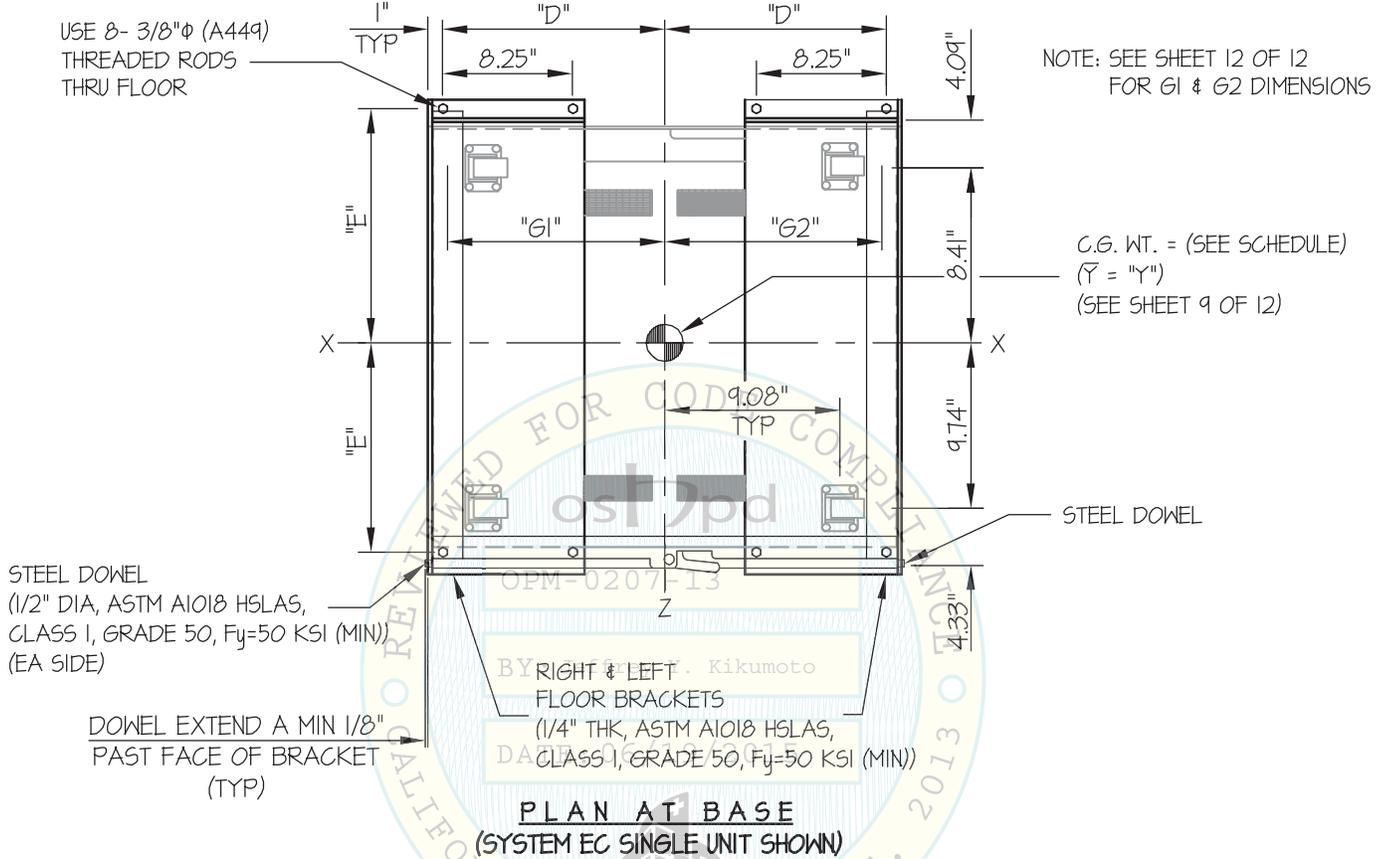

*Jonathan Roberson*

REGISTERED PROFESSIONAL ENGINEER  
JONATHAN ROBERSON  
No. 4197  
EXP. 6-30-2016  
6/18/15  
STRUCTURAL  
STATE OF CALIFORNIA

## CAREFUSION

## SINGLE STORAGE CABINETS

DES. **J. ROBERSON**

JOB NO. **11-1454**

DATE **6/18/15**

SHEET

**9**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE SLAB ON METAL DECK



MODEL	WEIGHT (LB)	"A" (in)	"Y" (in)	"B" (in)	"D" (in)	"E" (in)	* TU (lb)	* VU (lb)
CONSOLE UNIT	360	48.50	24.25	30.00	14.375	11.75	780	221
MED STATION	1050	46.75	23.38	22.75	10.75	13.25	2469	724
PYXIS STATION SINGLE UNIT	1460	78.06	39.03	30.30	14.53	13.25	4547	1020
SYSTEM 30 SINGLE UNIT	1460	78.06	39.03	30.30	14.53	13.25	4547	1020
SYSTEM EC SINGLE UNIT	1460	78.06	39.03	30.30	14.53	13.25	4547	1020

\* VALUES DO NOT INCLUDE  $\Omega$ .

NOTE: VALUES ARE FOR  $z/h \leq 1$ .

BY: Jeffrey Y. Kikumoto

DATE: 06/19/2015

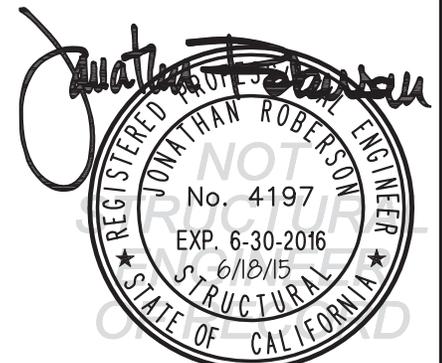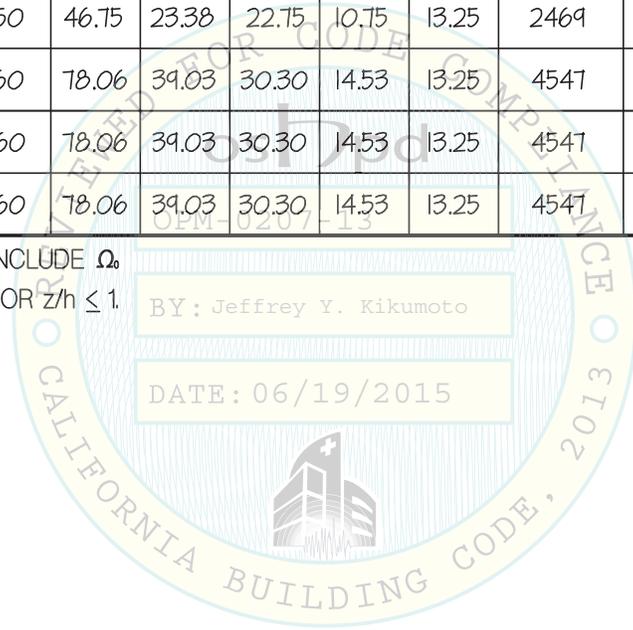

## CAREFUSION

DES. **J. ROBERSON**

SHEET

# 10

## SINGLE STORAGE CABINETS


JOB NO. **11-1454**

DATE **6/18/15**

OF **12** SHEETS

SEISMIC SUPPORTS & ATTACHMENTS

CONCRETE DETAIL

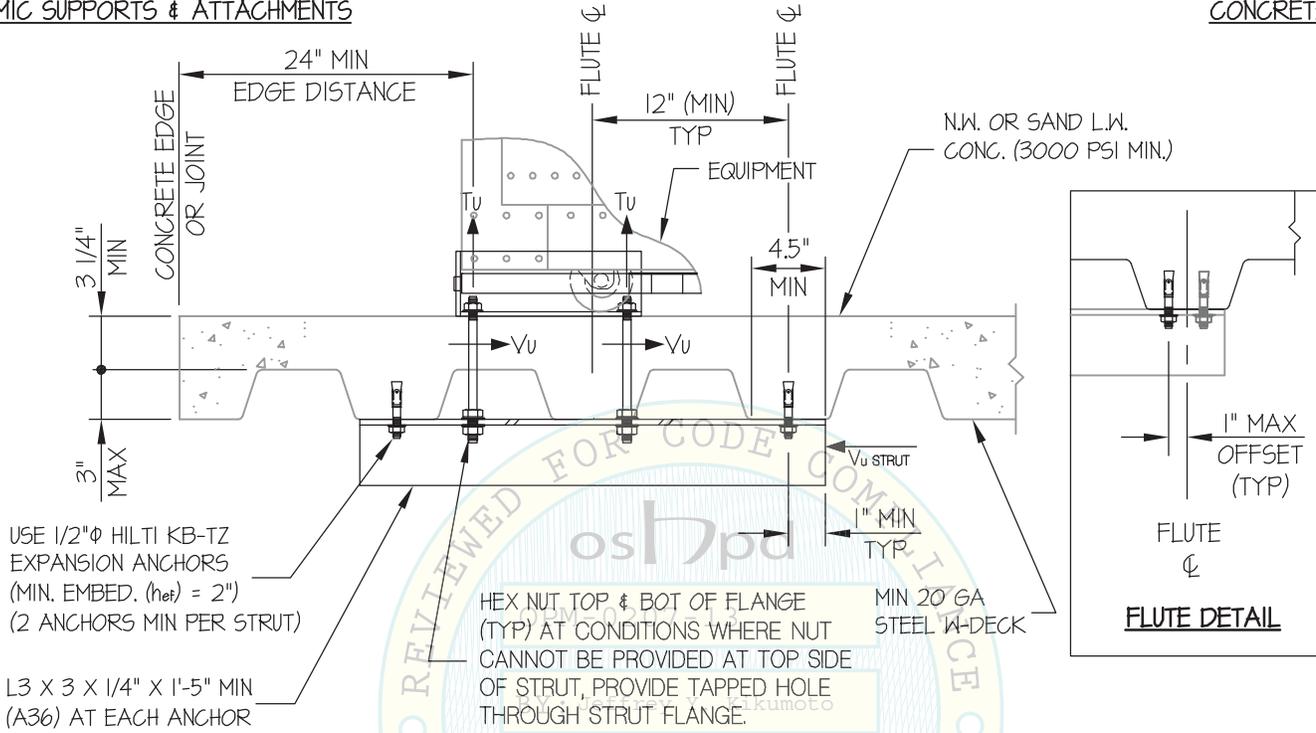

**MIN STEEL DECK REQUIREMENTS AND STRUT DETAIL**

NOTE: THERE IS A POSSIBILITY THAT 4  
THREADED RODS CAN BE ON A SINGLE STRUT.  
AN EXAMPLE CAN BE SEE ON PG 7 OF 12



## CAREFUSION

DES. **J. ROBERSON**

SHEET

**11**

## SINGLE STORAGE CABINETS

JOB NO. **11-1454**


DATE **6/18/15**

OF **12** SHEETS


SEISMIC SUPPORTS & ATTACHMENTS

FLOOR BRACKETS  
(1/4" THK, ASTM A1018 HSLAS,  
CLASS 1, GRADE 50, F<sub>y</sub>=50 KSI (MIN))


BASE DETAIL

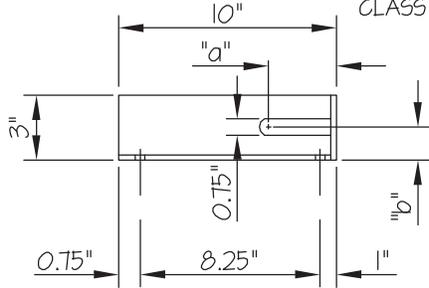

**SECTION A-A**

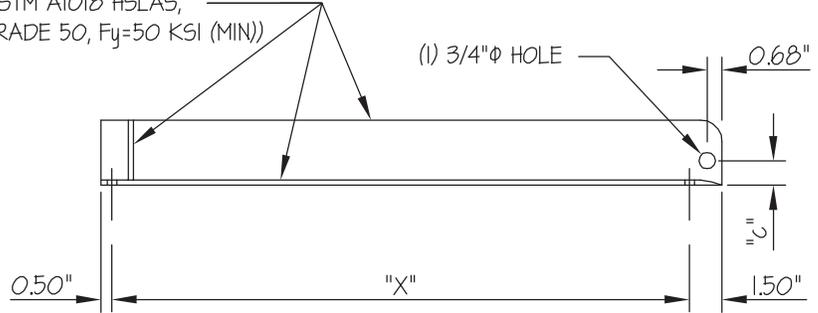


**SECTION B-B**

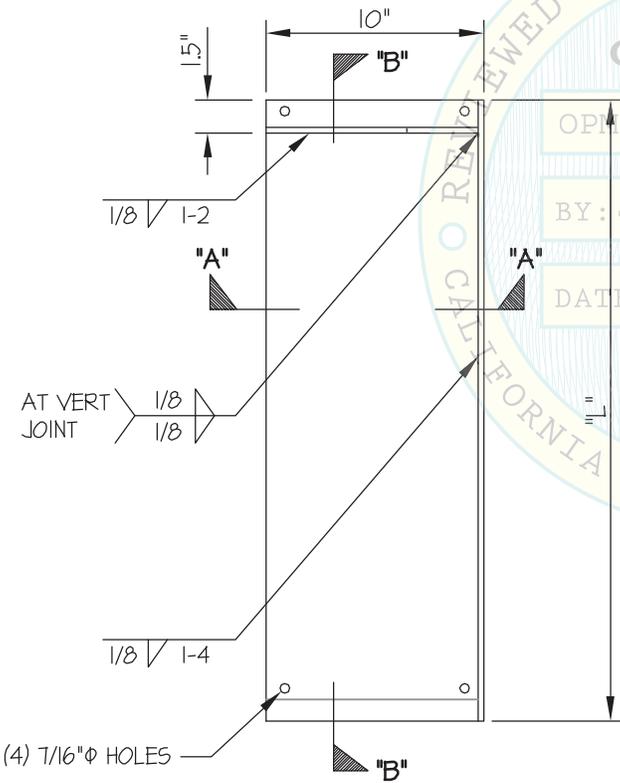

**PLAN**

(RIGHT FLOOR BRACKET SHOWN)

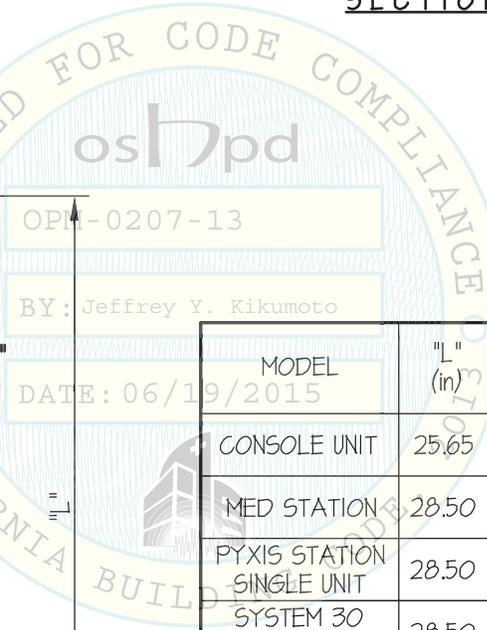

MODEL	"L" (in)	"X" (in)	"a" (in)	"b" (in)	"c" (in)
CONSOLE UNIT	25.65	23.65	3.12	1.795	1.55
MED STATION	28.50	26.50	3.12	1.425	1.18
PYXIS STATION SINGLE UNIT	28.50	26.50	3.50	1.885	1.61
SYSTEM 30 SINGLE UNIT	28.50	26.50	3.50	1.885	1.61
SYSTEM EC SINGLE UNIT	28.50	26.50	3.50	1.885	1.61

NOTE: THE SEISMIC ANCHORING KIT CONSISTS OF TWO FLOOR BRACKETS, A RIGHT FLOOR BRACKET AND A LEFT FLOOR BRACKET (SHOWN HERE) AND FRAME ASSEMBLY (SEE SHEET 12 OF 12) (TYP)



## CAREFUSION

DES. **J. ROBERSON**

SHEET

**12**

## SINGLE STORAGE CABINETS

JOB NO. **11-1454**


DATE **6/18/15**

OF **12** SHEETS

**SEISMIC SUPPORTS & ATTACHMENTS**

**BASE DETAIL**

STEEL DOWEL  
(1/2" DIA, ASTM A1018 HSLAS,  
CLASS I, GRADE 50, Fy=50 KSI (MIN))  
(EA SIDE)

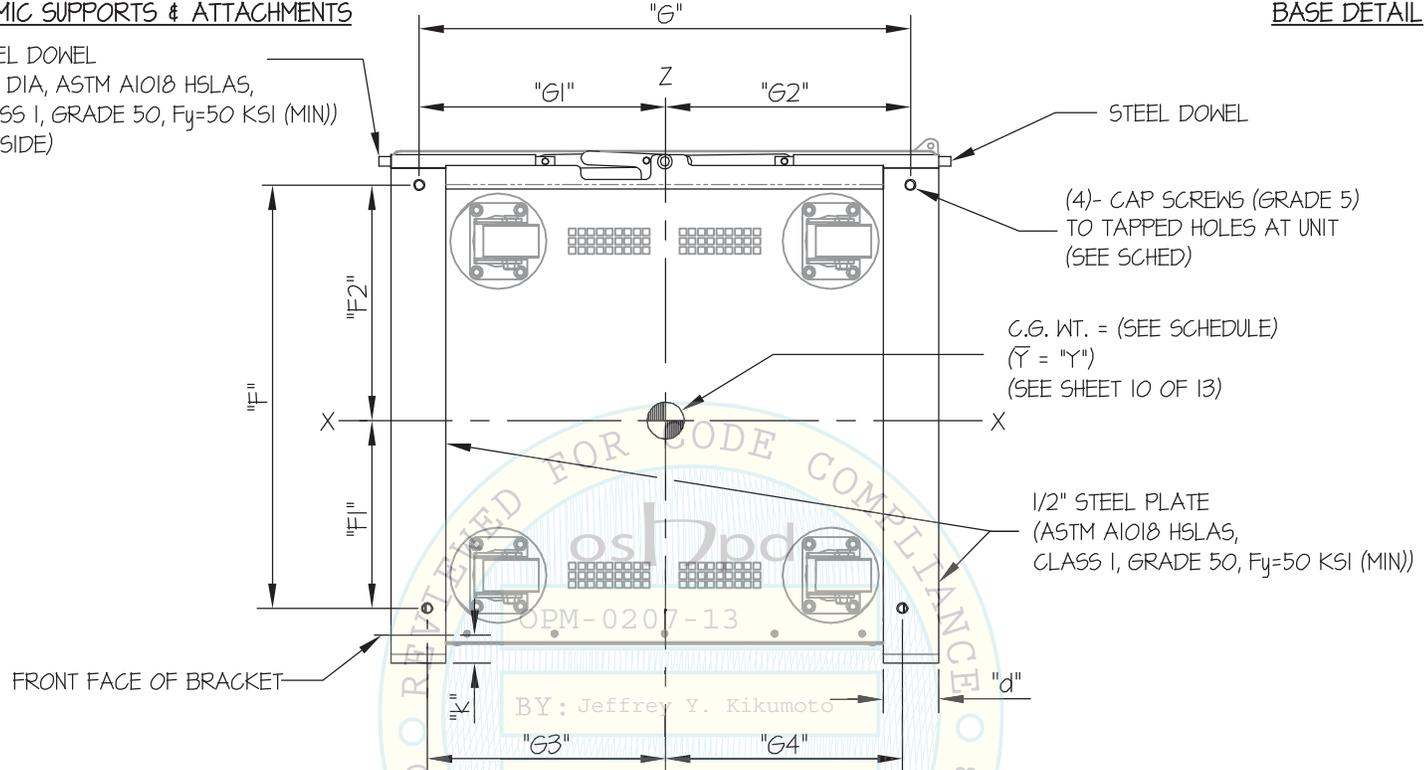

**REFLECTED PLAN AT FRAME ASSEMBLY**  
(FRAME ASSEMBLY SHOP ATTACHED TO UNDERSIDE OF UNIT BY MFG)

MODEL	CAP SCREW DIA (in)	"F" (in)	"F1" (in)	"F2" (in)	"G" (in)	"G1" (in)	"G2" (in)	"d" (in)	"G3" (in)	"G4" (in)	"k" (in)
CONSOLE UNIT	1/4	21.26	10.08	11.18	28.22	13.82	14.4	2.25	14.15	14.73	0.61
MED STATION	3/8	22.25	10.45	11.8	19.42	9.45	9.97	2.25	9.45	9.97	0.68
PYXIS STATION SINGLE UNIT	1/2	22.98	10.20	12.78	26.66	13.38	13.28	3.00	12.94	12.84	0.68
SYSTEM 30 SINGLE UNIT	1/2	22.98	10.20	12.78	26.66	13.38	13.28	3.00	12.94	12.84	0.68
SYSTEM EC SINGLE UNIT	1/2	22.98	10.20	12.78	26.66	13.38	13.28	3.00	12.94	12.84	0.68

**NOTE:** THE SEISMIC ANCHORING KIT CONSISTS OF TWO FLOOR BRACKETS, A RIGHT FLOOR BRACKET AND A LEFT FLOOR BRACKET (SEE SHEET 11 OF 12) AND FRAME ASSEMBLY (SHOWN ABOVE) (TYP)

